



EN

TS202*

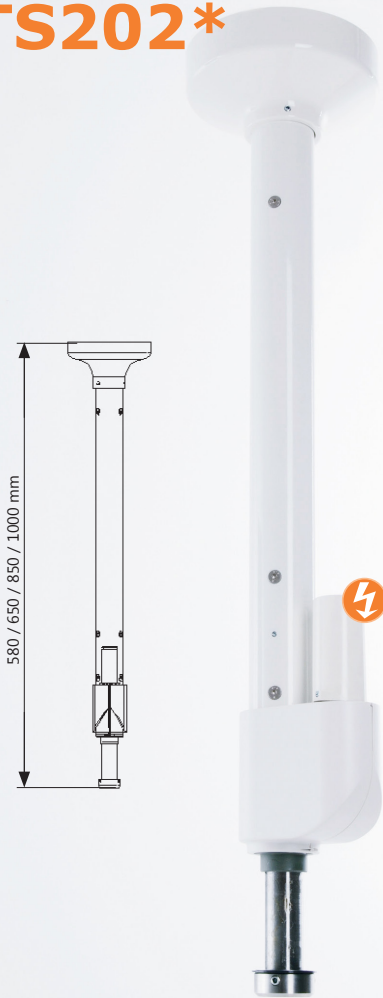

TS200*

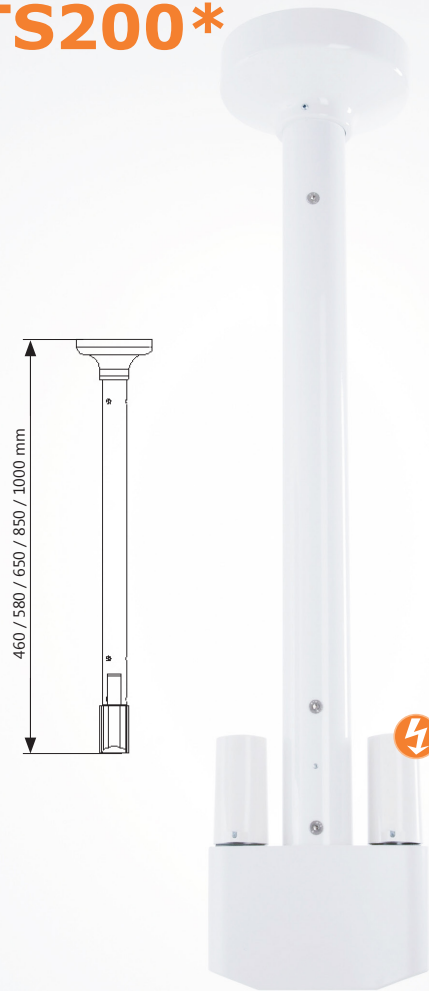

MAVIG

PORTEGRA²
360° Stationary Ceiling Column

PORTEGRA²
Twin Stationary Ceiling Column



Models: TS2020 - 23 Stationary 360° ceiling column (for column lengths refer to table below)
TS2001 - 05 Stationary Twin ceiling column (for column lengths refer to table below)

Made in Germany

IEC 60601-1

MAVIG Stationary Ceiling Columns – 360° and Twin

Many years of development has led to MAVIG's unique steel columns. The proven construction, sturdy design, and numerous tests ensure the highest level of safety available. The patented ceiling anchor provides an extra level of protection, stopping the system from falling in the event of a collision.

MAVIG columns are equipped standard with one electrical pin ⚡ for a direct connection to a lamp. If a solution with two electrical pins is necessary, the Twin column can be upgraded upon request.

- Each pin offers a 235° rotation (360° for the lower pin of the 360° column)
- One electrical connection: One 3-pole sliding contact ⚡
- Net load at spring arm adapter: Max. 18.0 kg (39.7 lbs) per pin (depending on length)
- Different lengths to match room height

Article-No.	360° Length*	Article-No.	Twin Length*
TS2020	1000 mm	TS2001	1000 mm
TS2021	850 mm	TS2002	850 mm
TS2022	650 mm	TS2003	580 mm
TS2023	580 mm	TS2004	650 mm
		TS2005	460 mm

* Ceiling column lengths are measured from the top of the mounting flange to the bottom of the column structure.

MAVIG GmbH
X-Ray Protection and Medical
Suspension Systems
Stahlgruberring 5
81829 Munich
Germany
PO Box 82 03 62
81803 Munich
Germany
Phone +49 (0) 89 / 420 96 - 0
Fax: + 49 (0) 89 / 420 96 - 200
e-Mail info@mavig.com
www.mavig.com
ТОВ „ОНКОМЕДІКА“
Уповноважений
представник в Україні
01014, Україна, м. Київ, вул.
Звіринецька, 63
+38 (044) 228 84 17
+38 (097) 959 10 89
info@oncomedica.ua
www.oncomedica.ua
www.mavig.com.ua