

Synergy RA632

ALL-ROUND PROTECTION: COAT

The simple, convenient way of putting it on adds to the comfort. The coat is simply closed with a hook-and-loop fastener at chest level and sits correctly right away.

Protection/Lead Equivalent Standard Collection

NovaLite, Leadfree
or Standard lead



Front:
0.25 mm Pb, 0.35 mm Pb
or 0.50 mm Pb

Back generally:
0.25 mm Pb

Protection/Lead Equivalent Levio Edition

NovaLite, Leadfree
or Standard lead



Total value Levio-Panel:
0.42 mm Pb or 0.50 mm Pb
(outside 0.25 mm Pb)

Back generally:
0.25 mm Pb

- ✔ Completely overlapping front parts achieve the selected lead equivalent value in the front
- ✔ Ideal fit, maximum comfort for the wearer
- ✔ Continuous support thanks to the wide elastic stretch belt (included)
- ✔ Standard equipment includes shoulder pads
- ✔ In the bottom area, the coat is additionally secured with a snap closure
- ✔ Type testing and certification in accordance with (EU) 2016/425, IEC 61331-1:2014 / DIN EN 61331-1:2016 and IEC 61331-3:2014 / DIN EN 61331-3:2016

For additional shoulder relief, the coat can be combined with our **BRYGGA** shoulder relief system.



Sizes*, recommended for:

XS	Chest girth ≤ 84 cm	Waist/hip girth ≤ 98 cm
S	Chest girth approx. 85 - 95 cm	Waist/hip girth approx. 99 - 108 cm
M	Chest girth approx. 96 - 106 cm	Waist/hip girth approx. 109 - 119 cm
L	Chest girth approx. 107 - 117 cm	Waist/hip girth approx. 120 - 130 cm
XL	Chest girth approx. 118 - 128 cm	Waist/hip girth approx. 131 - 141 cm
XXL	Chest girth approx. 129 - 139 cm	Waist/hip girth approx. 142 - 152 cm
3XL	Chest girth approx. 140 - 150 cm	Waist/hip girth approx. 153 - 163 cm


* The size depends on the girth of the chest, waist, or hip, depending on which of the three measurements is the largest.



A Safe Fit

Colour-coded (Titan) contrasting safety zones give you the certainty of being optimally protected in all situations.

If the front panel is within this clearly recognisable zone when the protective clothing is closed, then you are covering the full 60% of the girth with the maximum protective value.

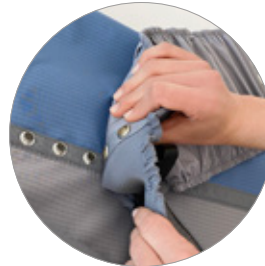

Snap Button Fastening

Optimal adjustment to your waist height is ensured by a variable snap button fastener of the belt. This allows it to be attached variably in height and to be buttoned to the coat on one or both sides.

Levio®

Optimised Protection Concept



BRYGGA Shoulder Relief



Outlast Temperature Regulation

Higher Protection – Levio Edition

The area with the reinforced protective material is clearly recognisable by its colour.

Safety Zone

This zone (colour-contrasted in Titan) ensures a risk free, seamless and safe closing of the coat RA632.

Back Relief

The wide elastic stretch belt makes it easier for you to keep your posture in daily routines.

Free Movement

For greater freedom of movement, the coat has a cut-out in the back. Of course, the knees are still optimally protected at the front and sides.



Outer Material

Colours
(see front page)



Optional:
ComforTex® HPMF,
ComforTex® HPMF Hybrid


Sizes

from XS to 3XL
(see front page)



Coat length / Body height:
90 cm / 150 - 160 cm
100 cm / 160 - 170 cm
110 cm / 170 - 180 cm
120 cm / 175 - 190 cm
130 cm / 185 - 195 cm
140 cm / 195 - 205 cm

Optional

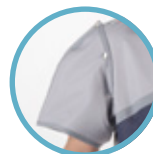

On request:
Outlast
BRYGGA
Sleeves
Breast pocket
Embroidery
Made-to-measure protection

Accessories/Replacement

Bolero
Thyroid/sternum protection



Stretch belt:
(Adjustment via
snap fastener)